

# McKESSON

## Durable Rollator

Provides a total weight bearing capacity of 300 lbs (136 kg).

For use with MFR#'s 146-10257BL-1 and 146-10257RD-1.

### ASSEMBLY INSTRUCTIONS

Do not attempt to adjust or operate Rollator without carefully reading all instructions in this manual.

1. First remove all contents from carton: backrest, handle with tube, Rollator frame and basket.
2. Open frame and slide padded backrest into frame brackets making sure snap bottoms protrude through holes.
3. Loosen tightening knob on handle receptacle and insert handle tube into frame. The brake cable should be positioned on outside of frame. Set handle to desired height. Tighten knob by turning clockwise until handle is securely in place. Make sure both handles are set to the same height before use.



### SAFETY PRECAUTIONS

1. Please make sure the Rollator is fully opened, seat is down and the folding bar under seat is locked in the down position.
2. Be sure that both knobs are fully tightened and are at equal heights. Test Rollator for stability before applying full weight to hand grips.
3. Test the loop locks before operating the Rollator. When wheels are locked (see illustrations for locking and releasing locks below), lock bar on rear wheels should prevent wheels from turning. If the wheels turn while brakes are locked, please take unit back to the dealer for adjustment. When locks are released, the wheels should turn freely.
4. Always lock loop locks before sitting and do not use seat when unit is on an incline or uneven ground.
5. Replace tires when significantly worn or damaged.


### IMPORTANT SAFETY NOTICE

Do not have anyone push you while you are seated on

the Rollator. This is a walking aid only and is not to be used as a transportation device.

### LOOP LOCK OPERATION

Always engage and disengage wheel locks on both sides simultaneously.

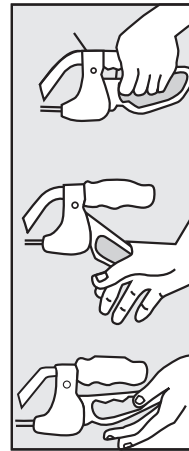

1. To Slow Down or Come to a Temporary Stop — Pull up slowly on the loop locks with four fingers until unit slows down or comes to a temporary stop.
2. To Temporarily Lock Wheels in Position — Press down on both loop locks with palms until you hear a “click” that will tell you wheels are locked.
3. To Release Wheel Locks — Push up on both loop locks with palms of hand to original position.

### INSTRUCTIONS TO ADJUST THE HAND BRAKE SYSTEM


#### IMPORTANT SAFETY NOTICE

The hand brake system fitted on this rollator complies with the manufacturer's quality standards. The hand brake system has been preset (adjusted) to the best performance settings. To better suit you personally, please follow the instructions set out below. You may need to adjust the hand brake systems from time to time as they will wear out with normal use.

Hold firmly in place do not turn.

To tighten the hand brake adjustment. “turn the nut clockwise”

To loosen the hand brake adjustment. “turn the nut counter clockwise”

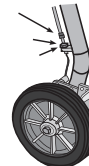

**WARNING:** This product can expose you to bisphenol A (BPA), which is known to the State of California to cause birth defects or other reproductive harm. For more information, go to [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov).

Technical support, replacement parts and warranty information? Call 1-866-694-5085 General product questions? Call 1-800-777-4908

#### Satisfaction Guaranteed

For complete details, please visit [mms.mckesson.com/mckesson-brands](http://mms.mckesson.com/mckesson-brands).

Distributed by McKesson Medical-Surgical Inc.  
Richmond, VA 23233  
PVN B0518  
Made in China